



Instrument Engineering Change Notice

NO. TEC6-2006

Sakura Finetek Japan Co. Ltd.

PRODUCT NAME : Tissue-Tek TEC 6 EM DATE : October 27, 2020
 PRODUCT CODES : 5104, 5108, 5112, 5128 ORIGINATOR : Yoshikawa
 IECN TYPE : SOFTWARE CHANGE REFERENCE NUMBER : 16877
 SERIAL NUMBERS AFFECTED : See Remarks ATTACHMENTS : TF10-20-1459-01_TEC6
 FIELD RETROFIT REQUIRED : YES Software update
 procedure manual.pdf,
 TF10-20-1461-01_Repair
 work report and
 completion
 confirmation.pdf,
 TF10-20-1470-00_Repair
 work report and
 completion confirmation
 for stock item.pdf,
 A-SAL2-0800-04-A.zip

Title: Software Update

[REASON AND DESCRIPTION FOR CHANGE]

To fix the incorrect date/time processing that occurs in March in a leap year, such as that the date calculation decreases by one day and that the invalid clock error is displayed, the control board program was modified.

With the change, the control program installation packages including the modified program A-SAL2-0100-04-A were updated.

[PART CHANGE] (IVDP=Inventory Disposal, ICHB=Inter Changeable)

Description	Old	QTY	IVDP	New	QTY	ICHB
Installation Package (5104, 5108, 5112, 5128)	A-SAL2-0800-03-A	1	YES	A-SAL2-0800-04-A	1	YES

[REMARKS]

J0: 5104 0197 and after
 A1: 5108 0334 and after
 E2: 5112 0244 and after
 J2: 5128 0101 and after

* For 5116 (JC2), IECN #TEC6-2007 will be issued separately.

[DISTRIBUTION]

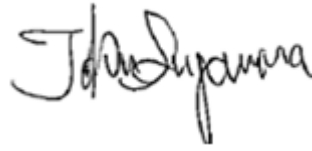
[APPROVAL]

SFA

Name : T. Sugawara

SFJ

Signature :

A handwritten signature in black ink, appearing to read 'T. Sugawara', enclosed within a thin vertical rectangular border.

SFE

SCN